



# OsmoTECH<sup>™</sup> Single-Sample Micro-Osmometer

Critical during development and manufacture of biological therapeutics and injectable solutions

### Supports 21 CFR Part 11 compliance

Osmolality testing is a robust in-process control and quality check you can implement throughout process development to ensure optimum media preparation, cell health, and product quality and yield.

OsmoTECH's data management features **support 21 CFR Part 11 compliance**, and its gold-standard testing method uses freezing point depression to provide the accuracy and precision you need to optimize your process and quality control program.



ADVANCED

# Why test osmolality?

therapeutic molecules, process development. reduce the risk of expensive batch to patients. — If you're not routinely testing osmolality, you're putting your biological drug at risk.

# Why OsmoTECH?

OsmoTECH is an easy-to-use, reliable system that provides a critical quality measure you can easily fit into your workflow.

#### With OsmoTECH, you can:

- · Optimize cell health and improve product yield
- Manage data and security with features that support 21 CFR Part 11 compliance
- Ensure accurate media and buffer preparation
- Identify process issues and avoid expensive batch loss
- Avoid product stability issues and patient discomfort

#### **Optimize your bioproduction process**

You need robust process development—where you can create reproducible processes at each stage-so you can have confidence you're developing a quality therapeutic. Testing osmolality gives you that confidence. At key points during processing, osmolality testing is critical-during upstream, downstream and as part of the quality check during formulation and fill.

# Your biotech production needs precision— **OsmoTECH** delivers

The bioproduction process is demanding. Timing is critical, and there are many inherent risk factors, such as human error, that impact your biological product development process. Each phase in the development process is critical and you need confidence that your therapeutic meets specification.


By the time your biological drug is ready for production release, you've invested time and money. You need confidence that your quality control is robust and your development process is reliable.

#### Why freezing point depression?

Freezing point depression is the gold standard osmolality testing method and the preferred standard throughout the world.

- Measurements are fast and accurate
- Easy-to-use and reliable
- · Works well with small sample sizes
- · Works with almost any liquid matrix
- Provides a full characterization of concentration

Other methods, including: specific gravity, refractive index, and conductivity, do not provide a complete view of a sample. Only freezing point osmolality is truly independent from the size, ionization status, shape, and other physical characteristics of the liquid solutions.

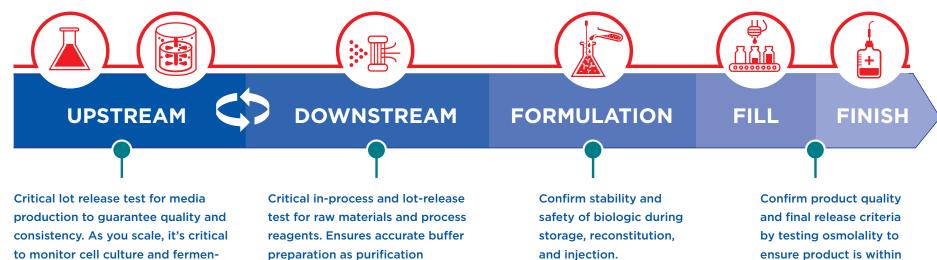

begins, and proves complete

buffer exchange to minimize

expensive batch failures and

lengthy quality investigations.

and yield.

tation to ensure optimal cell health


to optimize drug product quality

specific range.

### **Features and benefits**

#### Features that support 21 CFR Part 11 compliance

- Password-protection: Ensures security and traceability
- Audit trail: Time stamped records can't be modified
- Fraud detection: Completed test records can't be modified
- Access controls: Administrators can set up user accounts
- Electronic signatures: User passwords encrypted for data storage

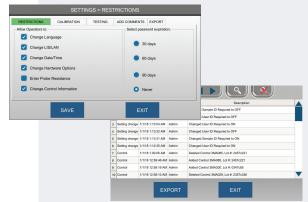

#### Manage and transfer data efficiently

- Storage: 1,000 test results and 10,000 events
- 2D barcode scanner: Eliminate transcription errors
- Download data: Supports .csv and .pdf formats: or disable
- Secure LIS link: Integrates with Laboratory Information System (LIS) using TCP/IP
- · Embedded web server: Remotely view, print, and download results from a network browser
- Network time synchronization: Ensure accurate time stamp/record alignment
- Dot matrix printer (optional): Prevents records from fading

#### Easy-to-use, flexible system

- Fast: 90-second test time
- One-step direct sampling: Simply aspirate sample and load
- Touchscreen: Color-coded, easy-to-use and intuitive
- Onscreen messaging: Instructional, status and statistics (mean, standard deviation, coefficient of variation)
- Clean room compatibility: Plastic corrugate consumable box
- Prevents sample transfer loss: Perfect for testing slightly viscous samples

## Request a quote today: aicompanies.com/OsmoTECH-quote

#### Parts and supplies

Part number	Product description
Instrument	
OsmoTECH	Single-Sample Micro-Osmometer
Reference Solutions	
3MA005	50 mOsm/kg calibration standard, 10x2 mL
3MA085	850 mOsm/kg calibration standard, 10x2 mL
3MA200	2000 mOsm/kg calibration standard, 10x2 mL
3MA029	Clinitrol™ 290 reference solution, 10x2 mL
3LA028	Osmolality linearity set, 5x2x5 mL
3LA011	100 mOsm calibration standard, 10x5 mL
3MA020	200 mOsm calibration standard, 10x2 mL
3MA040	400 mOsm calibration standard, 10x2 mL
3LA051	500 mOsm calibration standard, 10x5 mL
3LA091	900 mOsm calibration standard, 10x5 mL
3MA100	1000 mOsm calibration standard, 10x2 mL
3LA151	1500 mOsm calibration standard, 10x5 mL
Osmometer supplies and	d accessories
TECH500	Micro-Sample Test Kit: 2 boxes of 250 tips and cleaners; 1 plunger wire
3M0825	20 µL Ease-Eject™ Sampler
135022	Printer, dot matrix

Aqueous solution
20 ± 1 µL
90 seconds
Single sample
mOsm/kg H <sub>2</sub> O
1 mOsm/kg H <sub>2</sub> O
0 to 2000 mOsm/kg H <sub>2</sub> O
0 to 400 mOsm/kg $H_2O \le 2$ mOsm/kg $H_2O$ from nominal value (1 SD)
> 400 to < 1500 mOsm/kg $H_2O: \le 0.5\%$ from nominal value (1 SD)
≥ 1500 to 2000 mOsm/kg H2O: $\leq$ 1% from nominal value (1 SD)
0 to 400 mOsm/kg $H_2^{0}$ : SD $\leq$ 2 mOsm/kg $H_2^{0}$
> 400 to < 1500 mOsm/kg $\rm H_{2}O:$ coefficient of variation $\leq 0.5\%$
$\geq$ 1500 to 2000 mOsm/kg $\rm H_{2}O:$ coefficient of variation $\leq$ 1%
$<1mOsm/kg~H_2O$ per 5°C (9°F) ambient temperature change
USB 2.0 Type A ports (2), USB 2.0 Type B ports (1), Ethernet 10/100, RJ45 connector port (1), dot matrix printer (optional)
Simplified Chinese, Czech, Danish, English, French, German, Greek, Italian, Japanese, Korean, Portuguese, Russian, Slovak, Spanish, Swedish, Turkish
-20°C to +45°C (-4°F to +113°F)
100 to 240 VAC (50/60 Hz)
60 Watts
38 cm x 36 cm x 29 cm (15" x 14" x 11.5")
6.0 kg (13.3 lbs.)
11.4 kg (25 lbs.)
One-year limited warranty on workmanship and parts

1. Subject to change

**Specifications**<sup>1</sup>

 Accuracy and precision (within run) specifications apply to Advanced Instruments standards and reference solutions. Performance at Reference Conditions: 20°C to 25°C (68°F to 77°F); 40 to 60% relative humidity

 Operating Conditions: Temperature 18°C to 35°C (64°F to 95°F); 30 to 80% relative humidity (non-condensing)

4. Dimensions when Micro-Sample Test Kit is on the instrument

©2018 Advanced Instruments. Clinitrol, Ease-Eject, and OsmoTECH are trademarks of Advanced Instruments. All other trademarks are the property of their respective companies. Advanced Instruments certifies that the technical features needed for 21 CFR Part 11 compliance are built into OsmoTECH. It is your responsibility to implement the necessary controls in your laboratory to comply with 21 CFR Part 11 requirements.

MP00042, Rev.1

# ADVANCED INSTRUMENTS

RoHS

Not for patient diagnostic use

.

Two Technology Way /781-320-9000 Norwood, Massachusetts 02062, USA 800-225-4034 Fax:781-320-8181 aicompanies.com

# Request a quote today: aicompanies.com/OsmoTECH-quote

PCN00667

CE